



Android TRIO WAVELENGTH DIODE LASER



CONTENT

01

Android system

02

Treatment Head Advantages

03

Specifications

04

Before & After

05


About us

Android system interface

Diode laser 13.3&15.6-inch 4K Android, wifi, Bluetooth, 3D animation smart
ipad screen, you can even Google, watch TV, listen to music, and take
screenshots.

- ▲ Remote control system
- ▲ Clinical patient database
- ▲ Leasing business system

Remote control system

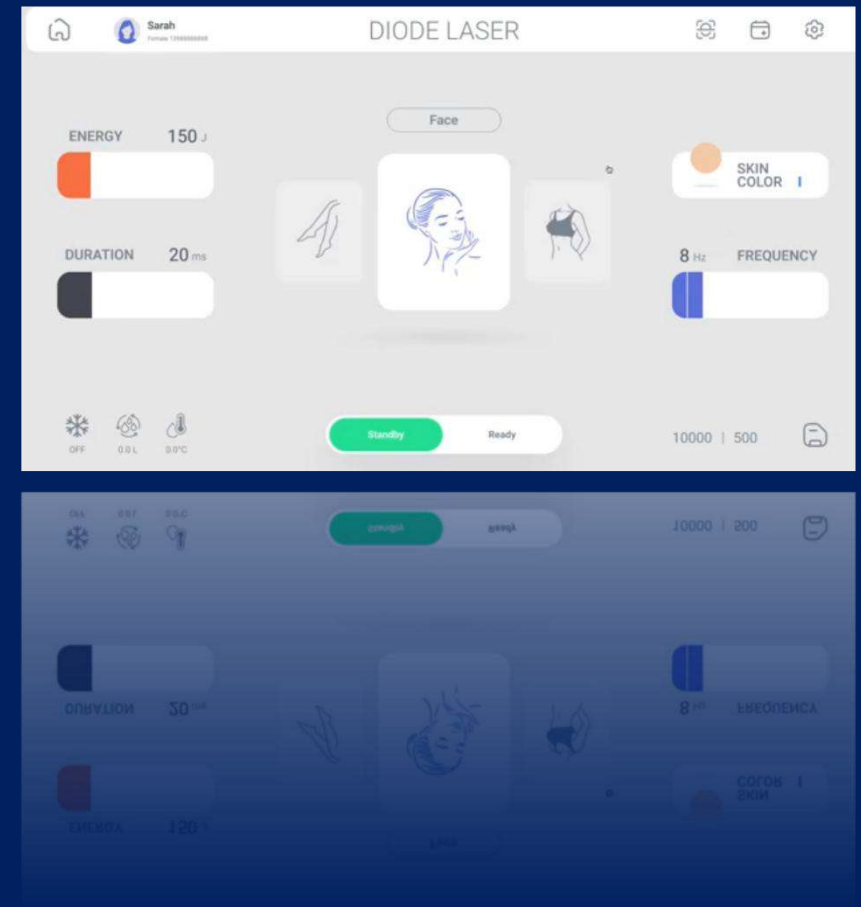



Only need a mobile phone

It is convenient for agents with leasing/installation business to remotely control the machine (Able to follow up the time or number of sending, turn on/off the machine remotely at any time, and understand the usage of the machine)

Remote control system



01

It is convenient for agents with leasing/installation business to remotely control the machine(Able to follow up the time or number of sending, turn on/off the machine remotely at any time, and understand the usage of the machine)

Remote control system Cloud data transfer
User Management system for patient database



04



02

Facilitate in-store management of large terminal beauty salons

--Information collection: Intelligent storage of usage parameters, customer number search, mass data retrieval
--Intelligent interface, set the total energy according to the area, convenient and quick to achieve hair removal effect




16 languages options

Customizable welcome page and screen logo



03

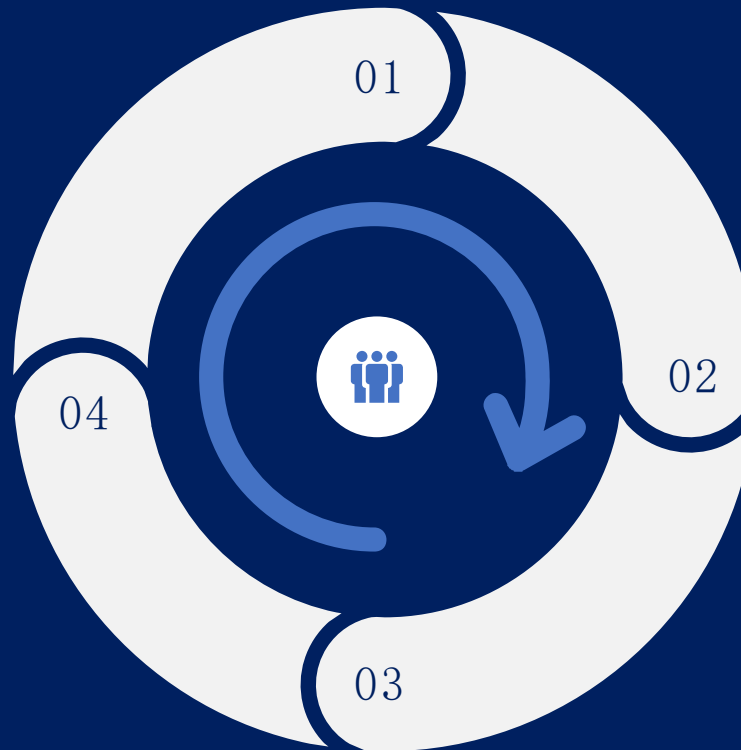

new HAIR DETECTION FUNCTION

File and save - data recording function




It is more conducive to recording customer cases and providing convenience for hair symptom analysis. Analyzing the skin and hair condition of the damage and nurture along with skin oil, normal, dry, dark and so on skin and hair types.



Theory

The principle of diode laser hair removal treatment is based on selective photothermal decomposition. The melanosomes in the hair follicle can selectively absorb the energy of the laser. The laser energy emitted by the machine is easily absorbed by the colored hair follicle without damaging the epidermal tissue. The energy emitted by the laser is absorbed by the hair and the pigment in the hair follicle.

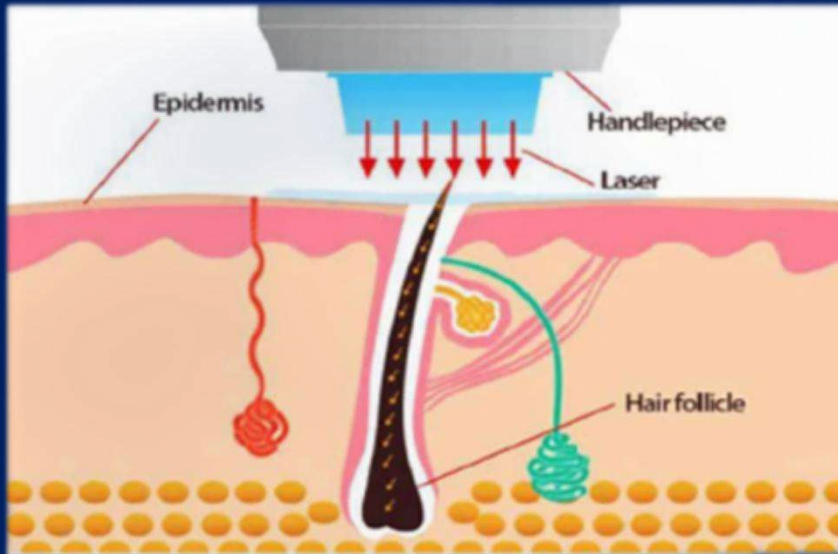

It is converted into heat, thus increasing the temperature of the hair follicle. When the temperature rises to a certain level, the hair follicle will be irreversibly damaged, and the hair will lose its original environment and be completely removed.



Straight hair

Curly hair

Fine hair

ALL HAIR TYPE APPLY



16mm
24mm

Big spot size

High power:

500w/800w/1000w/1200w/1600w/1800w/2000w

Handle shots 50 Million shots

Handle with 4k screen



755nm



808nm



1064nm



Suitable for all kinds of skin

755nm for white skin (fine, golden hair removal)

808nm for yellow neutral skin (brown hair removal)

1064nm for black skin (black hair removal)

Three WAVELENGTHS

755NM


For Light Skin

808NM

For Yellow Skin

1064NM

For Darker Skin

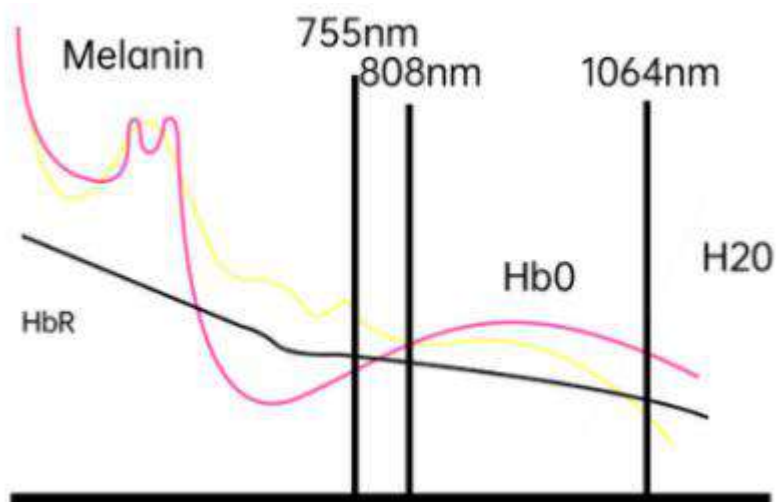

Suitable for all skin



FACIAL TIP

Improved reach easily treat areas that are normally hard to reach, including the ears, nostrils and glabella.




The LED screen displays healing energy and the number of rounds used, for easy viewing at any time



50 Million Shots

Guarantee 50 Million light emits under laboratory conditions. Under usage conditions, it can ensure effective light emission of 20 Million shots without energy attenuation

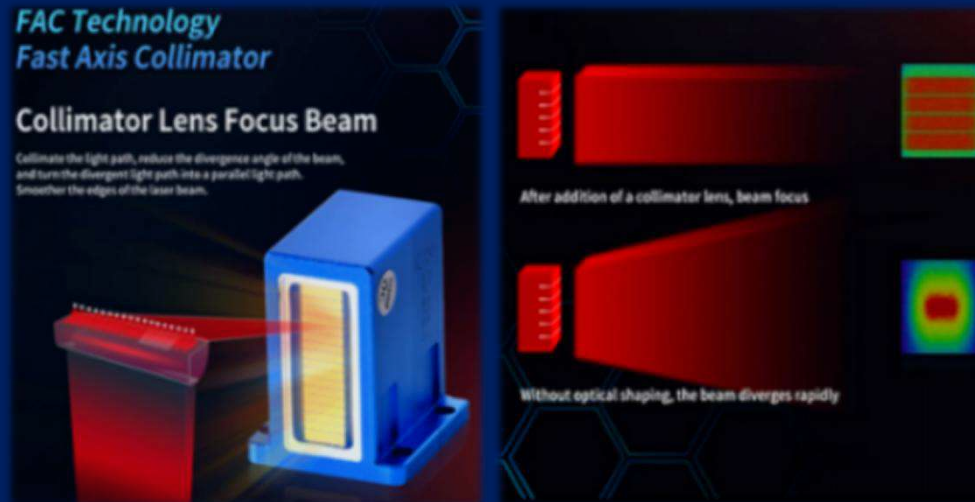

COOLING SYSTEM

ICE POINT PAINLESS

-16°




4K 13.3 & 15.6 INCH

LCD SCREEN

HIGH SENSITIVITY, HUMANIZED DESIGN, EASIER OPERATION.

Patient database

Viwe background data remotely

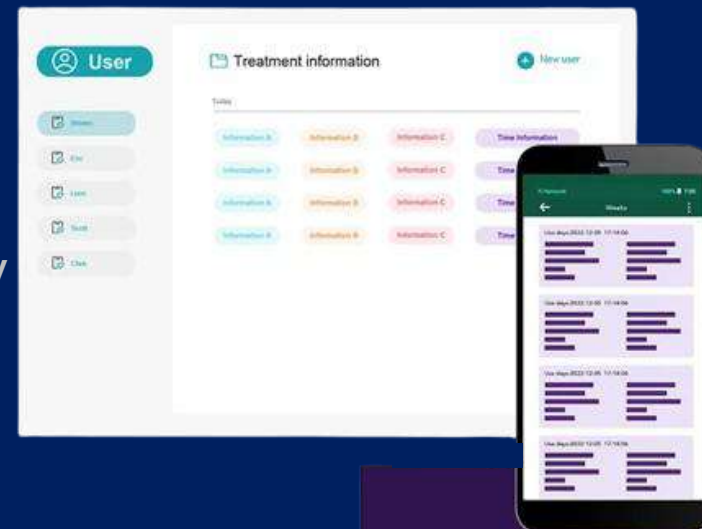


Before & After





755nm

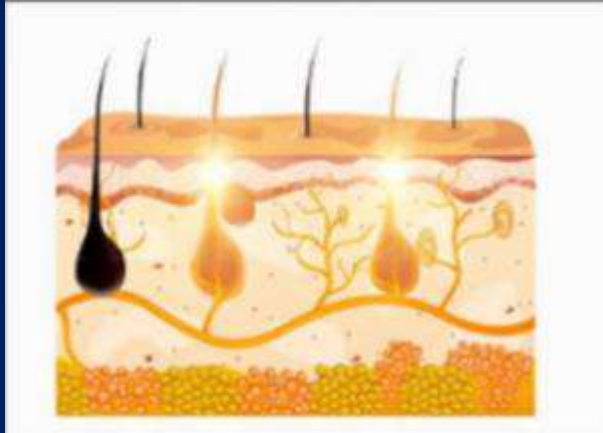


1-2 TIMES

Hair reduced by 75%
growth slows down



808nm

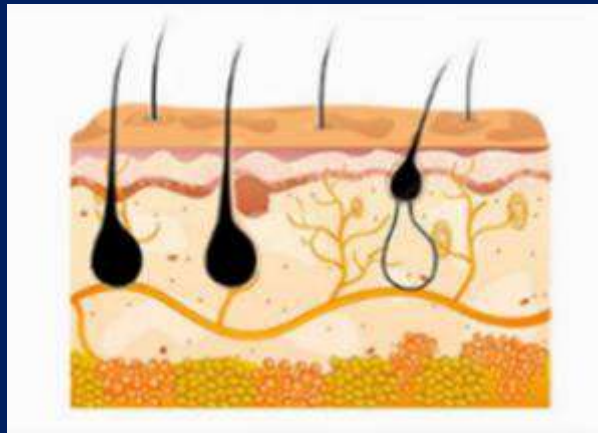

3-4 TIMES

less hair,
tender and smooth skin



4-5 TIMES ABOVE

Hair stops to grow as need



1064nm



Technical specification

Model Name	DepiMEDiX_Andriod
Machine Power	2000W
Laser Stack Power	500W
Diode Energy Density	1-150J/cm*
Diode Spot Size	16X24mm ² /φ6mm
Frequency	1-10Hz
Electrical Requirements	110/230VAC,15/20A max.,50/60Hz
Display	13.3&15.6"Android Screen
Net Weight	60KG
Dimensions	67*60*125cm
Stand-by Working	Continuously 24 hours



TECHNOLOGY OF BENEFIT

HIGH power TEC



ITALIAN water pump




MW power supply




The machine parts are selected from **international first-line brands** only to create high-quality machines.

AMERICAN ORIGINAL

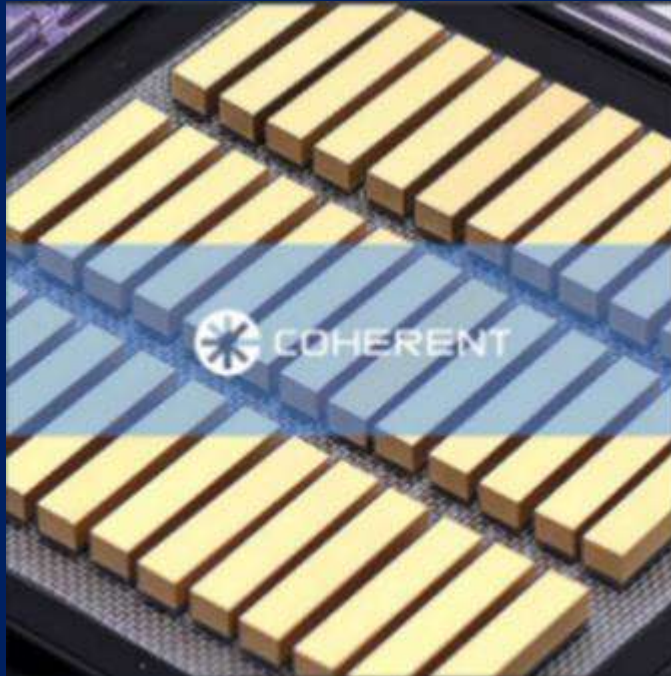

American Original Import

By Original American Coherent lasertechnology, the original bars are adopted in the laser generator to ensure that the output is stable and sufficient. The high quality bars guarantee a super long life span of the laser.

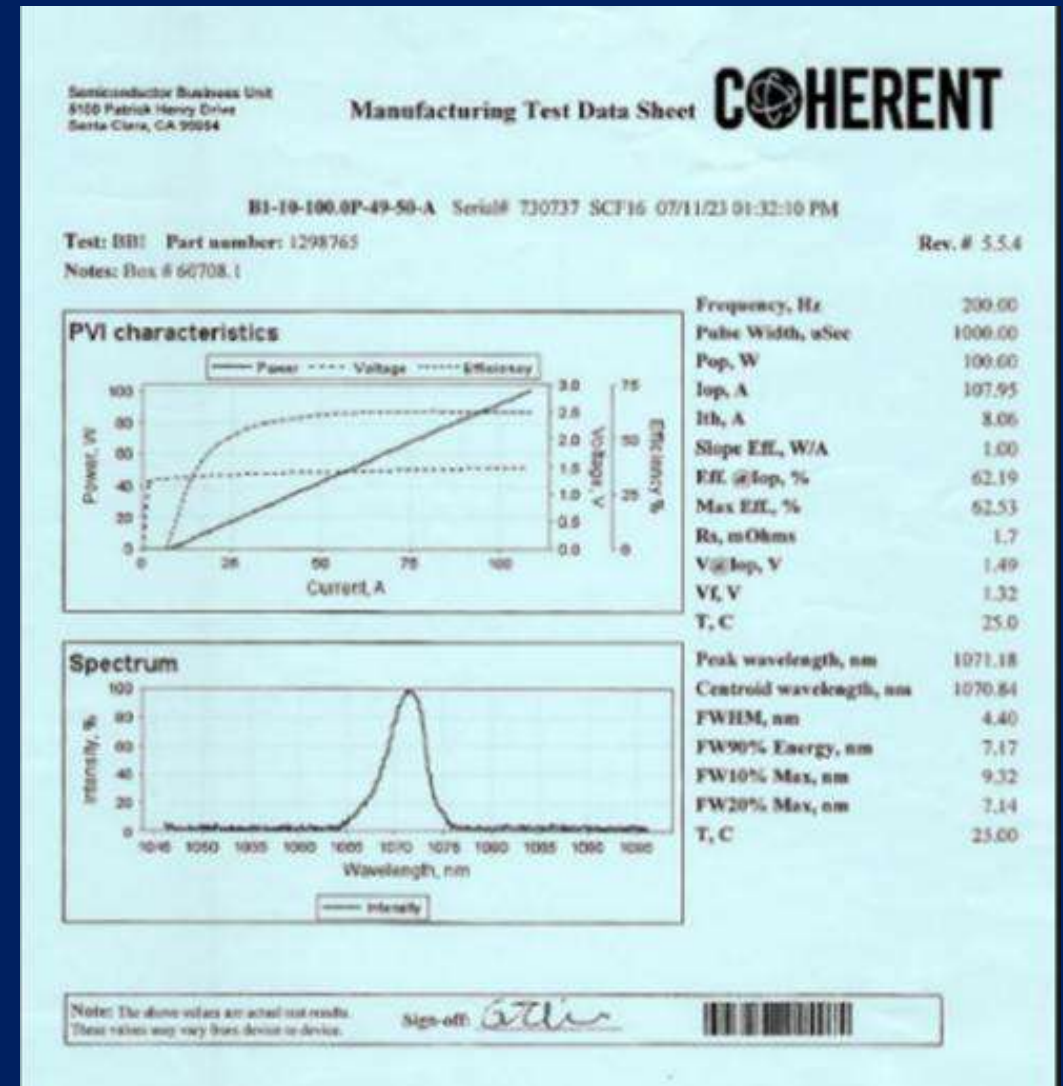

Views of TruMEDiO7





Welcome to visit

our Factory!




Our Certificates

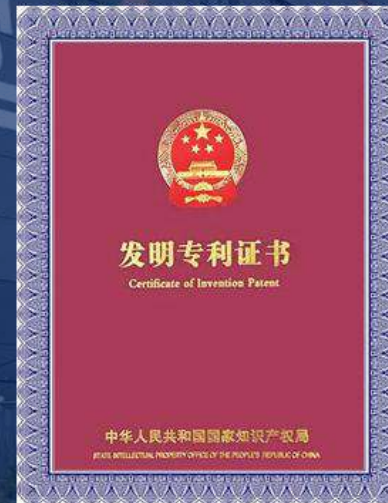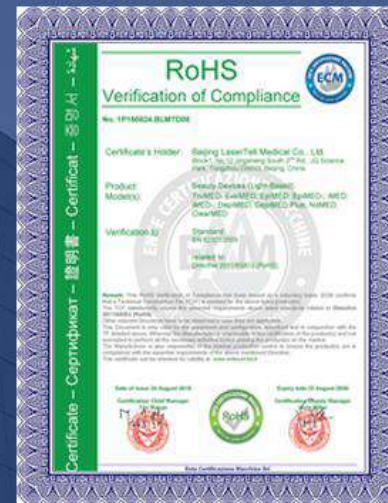

Well-recognized by Customers around 80+ Countries abroad with LaserTell Devices



Tailor-Made OEM/ODM Service

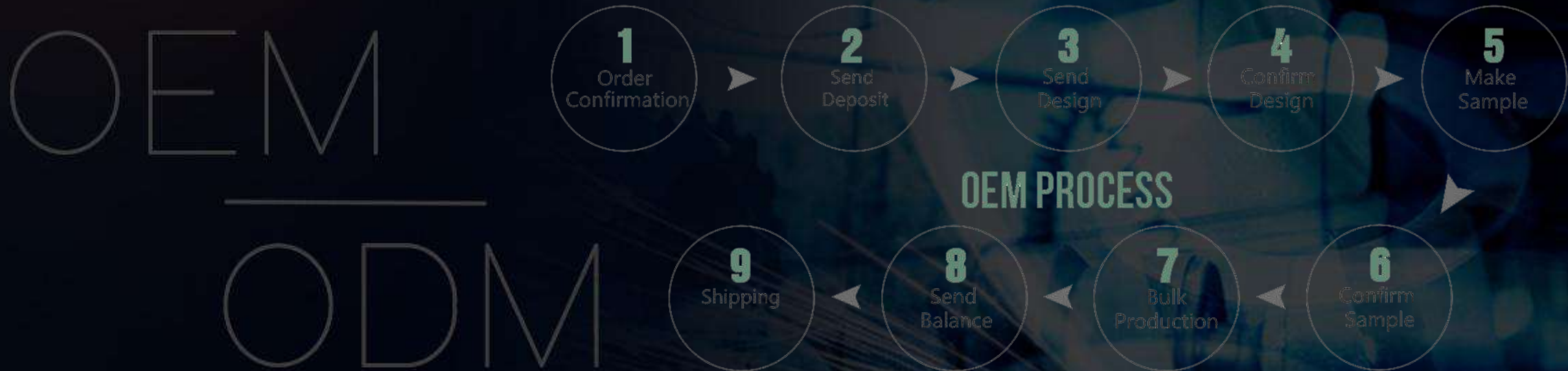


LaserTell has strong research and development team, including electronics engineers (designers, AVR programming; circuit layout) construction engineer, UI designer engineer QC control engineer...etc. We provide OEM/ODM design and program)

After-sale Service

2 years warranty available (2 + 1 charge) ex factory under normal and proper operation; free labor and material but without consumables parts.

Why 2+1 years warranty? (7 Simple Reasons)

1. Normally almost all the medical and aesthetic lasers with only 5 years lifespan, while the first 3 years' guarantee and full maintenance are right the key points.
LaserTell offers 2 years standard warranty and 1 year extension warranty. (Free-2, charge-1)
2. Lifelong Service (no matter 3 years or 5 years, or more...)
3. LaserTell, similar to a local service provider, has a full complement of certified technicians and engineers to assist you with all of your repair and maintenance needs, offering you faster service at much lower prices than the manufacturers
4. 24 hours + 7 days after-sale service and technical supports
5. **Service Manual** and Safety Instructions available
6. There is a comprehensive training last for 5-7 days in our factory and clinics, including how to operate, maintain and adjust the treatment parameters for different treatments, etc.

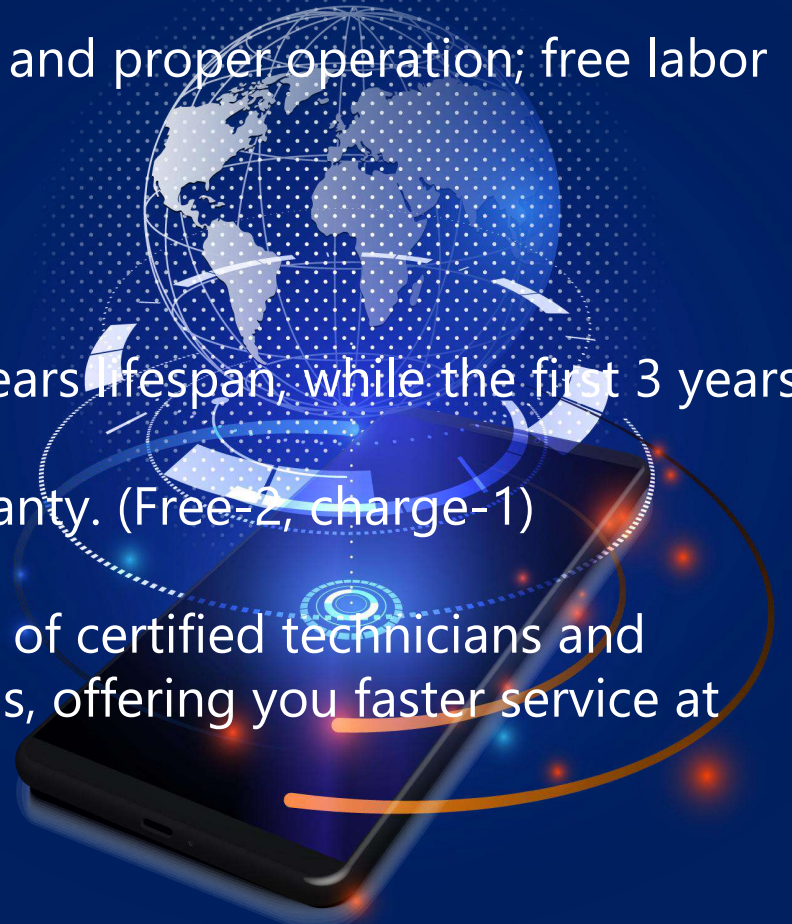

Thank you

