



#### CORE COMPETITIVENESS

# INDUSTRY LEADING TECHNOLOGY

OUR TECHNOLOGY IS NOT ONLY RECOGNIZED IN THE DOMESTIC MARKET, BUT ALSO ENJOYS A GOOD
REPUTATION IN THE INTERNATIONAL MARKET
CONTACT US














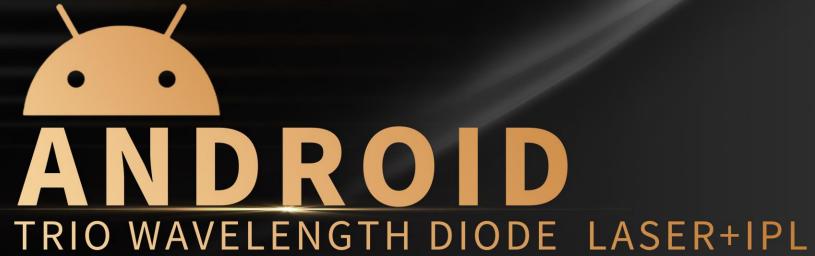




# Start Rental Business

- According to the lease agreement set the corresponding lease days and lighting times
- By setting the expired days of rental and the emitting shots to control whether the machine is running, more convenient for your leasing business.

**VISIT MORE** 

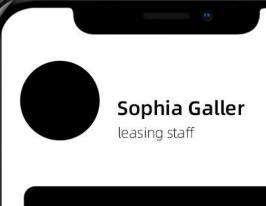

Remalning days:


2 9 9 9 9 9 2 8 8


Renral days

Number of ahots

# Patient database

- Remotely view User treatment data through mobile phones, and obtain effective information through data analysis.
- The patient's treatment parameters can be saved by the machine with one click and uploaded to the mobile patient database, and the patient data can be directly retrieved during the next treatment

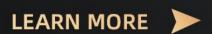

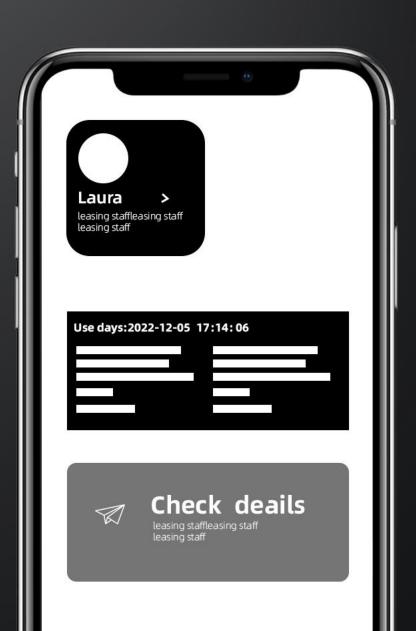

# **CUSTOMIZATION IN 16 LANGUAGES**

Built-in multiple languages to make the operation smoother more languages are being added...













NoW anytime, anywhere Only need a mobile phone to control the machine



super long warranty period Contains main components



**Remote Service** 



Maintenanceservice



**2YEARS WARRANTY** 























# SKIN DETECTION FUNCTION

Support skin detection function, real-time tracking treatment effec



New features, stay tuned...



200x lens

Scalp hair follicle amplification and fine detection

50x lens

Detect scalp skin condition

# 50 million shots


We always adhere to the use of real USA coherent lasers, enhance customer trust

Hot 800W/1200W/1600W/1800W/2000W

**USA** coherent laser

Learn more



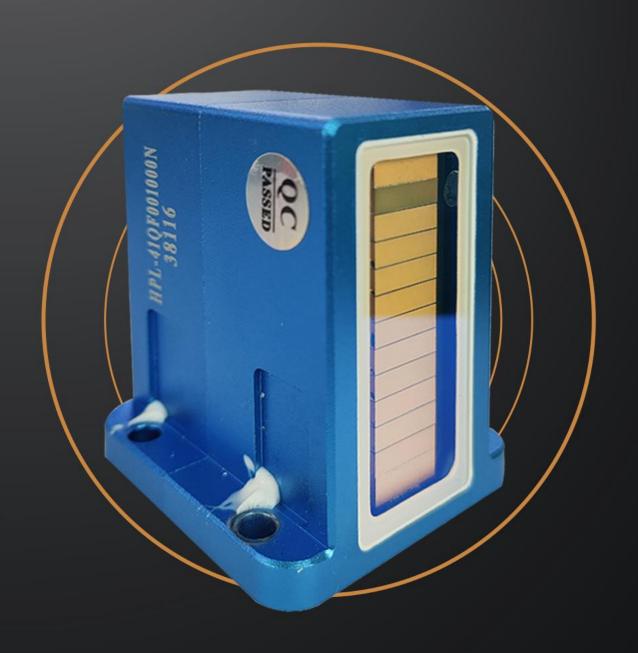





# Handle

**Multifunctional Detachable Design** 

# IPL HANDLE

#### Multifunctional beauty equipment

#### **OPTIONAL SPOT SIZE:**

8\*40mm<sup>2</sup>

#### **LMPORTED UK LAMP LONG SHOTS TIME:**

300,000--500,000 times

#### 6 FILTERS:

420-1200nm

430-1200nm

480-1200nm

530-1200nm

560-1200nm

640-1200nm



# 300,000-500,000 times

IPL handle long shots time



**Imported** 

#### **UK LAMP**

fo pi handlo af tho mschmo edoptsmpord Вгй sh lampa.whch can bousod for a long timo and omtt l ghtmany timas to ansuro a stabile woringstato.

Learn mons

# FILTERS

- 🕕 Permanent Hair Removal(HR)
- Skin Rejuvenation(SR)
- 🕕 Vascular Removal(VR)
- 🥨 Skin Tightening
- 🕕 Wrinkle Removal
- Breast Lifting
- Acne Removal
- Pigment Removal
- 🕛 🛮 Age Spot Removal etc.







STABILIZE








# HAIR REMOVAL PRINCIPLE


MAKE YOUR BEAUTY BUSINESS MORE EASIER



THE LASER SELECTIVELY ACTSON THE MELANIN IN THE HAIRFOLLI-CLE, DESTROYING THE-HAIR-GROWING AREA

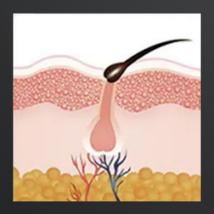

NATURAL HAIR LOSS, PREVEN-TREGENERATION, AND ACHIE-VETHE PURPOSE OFHAIR REMOVA

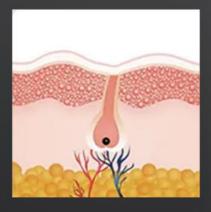

HAIR REMOVAL ALSO STIMULATES COLLAGEN REGEN-ERA-TION, SHRINKS PORES ANDMAKES SKIN FIRM ANDS-MOOTH

#### 755NM 808NM 1064NM WAVE ADVANTAGE

SPEED 808NM: GOLDEN STANDARD WAVELENGTH FOR ALL HAIR TYPES ALEX 755NM: SPECIFIC EFFECTIVE FOR BLONDE HAIR, FINE HAIR

LONG PULSE 1064NM:SPECIFIC EFFECTIVE FOR DARK, TANNED SKIN BEST QUALITY IPL SHR OPT LASER/ ND YAG /ALEXANDRITE LASER HAIR REMOVAL MACHINE PRICE






### MULTIPLE LANGUAGES AVAILABLE

### 4K 13.3 INCH

#### LCD SCREEN

HIGH SENSITIVIY, HURANIZED DESIGN. EASIER OPERATON.

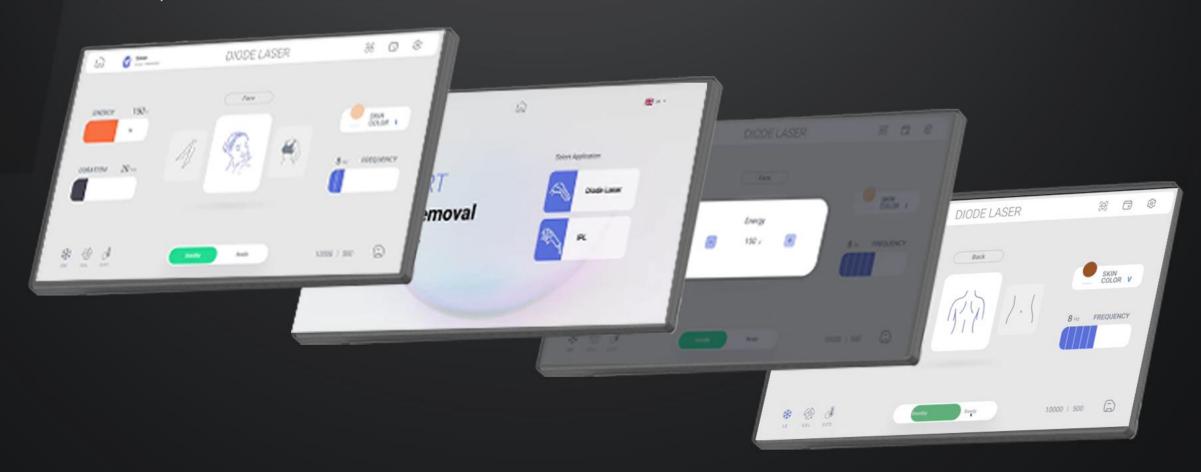

# TECHNICAL SPECIFICATION

MODEL NAME	TRUMEDI07 ANDRIOD 808+IPL
MACHINE POWER	4000W
LASER STACK POWER IPL POWER SUPPLY LPL SPOT SIZE	1000W 2000W 8X40MM <sup>2</sup>
DIODE ENERGYDENSITY  IPL ENERGY DENSITY  DIODE SPOT SIZE  IPL SPOT SIZE  IPL FILTER	1-150J/CM <sup>2</sup> 1-30J/1-50J 16X22MM <sup>2</sup> /Φ 6MM 8*40MM <sup>2</sup> 420-1200NM, 430-1200NM, 480-1200NM,
IFLFILIEK	530-1200NM, 560-1200NM, 640-1200NM

FREQUENCY 1-10HZ

ELECTRICAL REQUIREMENTS 110/230VAC,15/20A MAX,50/60HZ

DISPLAY 13.3"ANDROID SCREEN

NET WEIGHT 60KG

DIMENSIONS 67\*60\*125CM


STAND-BY WORKING CONTINUOUSLY 24 HOURS

# THREE WAVELENGTHS

755NM FOR LIGHT SKIN


808NM FOR YELLOW SKIN

1064NM FOR DARKER SKIN

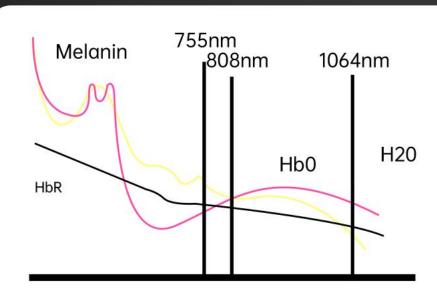



## AMERICAN ORIGINAL

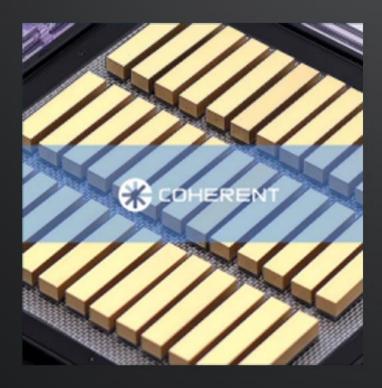

#### American Original Import

By Original American Coherent lasertechnology, the original bars are-adopted in the laser generator toensure that the output is stable and suffi-cient. The high quality barsguarantee a super zlong life span of the laser.

Semiconductor Business Unit 5100 Patrick Henry Drive Senta Clara, CA 95054




#### BI-10-100.0P-49-50-A Serial# 730737 SCF16 07/11/23 01:32:10 PM

Test: BB1 Part number: 1298765

Rev. # 5.5.4

Notes: Box # 60708.1

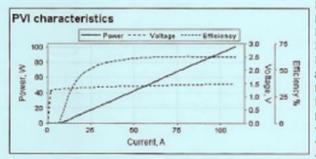

Frequency, Hz	200.00				
Pulse Width, uSec	1000.00				
Pop, W	100.00				
op, A	107.95				
th, A	8.06				
Slope Eff., W/A	1.00				
Eff. @Iop, %	62.19				
Max Eff., %	62.53				
Rs, mOhms	1.7				
V@Iop, V	1.49				
Vf, V	1.32				
r, c	25.0				

	100											11
Ŗ.	00						1	1				
ě	00						1	1				
928	40						1	1				
Intensity	20	-				1	/	1				
		Make	-		-	1		-	-	docuto	-	-
		1045 10	150	1055	1000	1000	1070	1076	1000	1005	1000	1005
						Wa	veleng	th, nm				

1,6	25.0
Peak wavelength, nm	1071.18
Centroid wavelength, nm	1070.84
FWHM, nm	4.40
FW90% Energy, nm	7.17
FW10% Max, nm	9.32
FW20% Max, nm	7.14
T, C	25.00

Note: The above values are actual test results. These values may vary from device to device.

Signoff: William



### **TECHNOLOGY OF BENEFIT**

ACTING ON THE SUPERFICIAL LAYER OF THE SKIN, COMBINED WITH 20HZ SLFHR® FLASHES, EFFICIENT, STABLE, PAINLESS.



#### WHAT IS THE ARESMIX®?

proprietary hair removal innovation is clinically proven method of laserhalr removal that is virtually painless, and it is the safe laser hair removal methodfor all skin types,including darkand tanned skin,lts unique qradualheating andrapid/fast pulse delivery method effectively damages the hair follicle and hindersre-growth, whilepreventing iniury to the surrounding tissue

# **UPGRADE TO4 WAVELENGTHS**

755NM

The 755Nm alexandrite performs hairremoval on Caucasian to olive skintones

808NM

The 808Nm gold standard inPermanent Hair Removal and issuitable on all pigmented hair andskin types

1064NM

The Yag 1064 wavelength is characterizecby lower melanin absorption, Specialized for darker skin types. making it a focused solution for darker skin types

# **TECHNOLOGY OF BENEFIT**

**HIGH** power TEC



**ITALIAN** water pump



**MW** power supply



The machine parts are selected from **international first-line brands** only to create hich-guality machines.



# Push data remotely



#### Energy

Piease erter a value from 1-300



#### **Frequency**

Please erter a value from 1-10



Online demo

**Leam more** 



Stable current management system

MW power

Learn about power





# Medical grade water flow precision


Italian water pump


Learn about water pump (







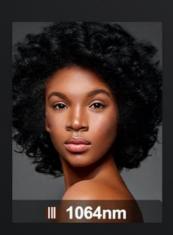


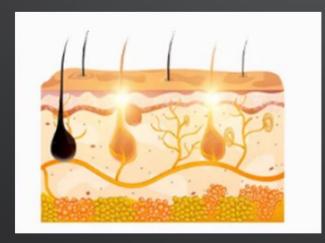

1-2 TIMES

Hair reduced by 75% growth sliws down



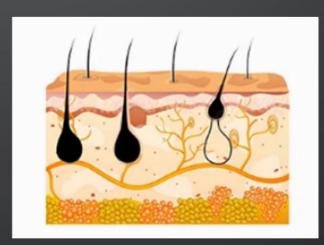

3-4 TIMES
less hair, tender and smooth skin





4 TIMES ABOVE
Hair stops to growthuse as need



