




HydroCool_PicoNEO



Picosecond

Tattoo Removal Machine

532nm, 755nm, 1320nm, 1064nm

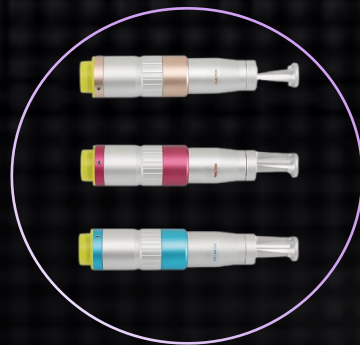

KOREA IMPORTED 7 LIGHT GUIDE ARM

360°

Inside placed lens mirror for better result



2-10mm adjustable
spot size




Weight balance hammer keeythe arm
always be safe



1064/532nm probe

4/532nm prob

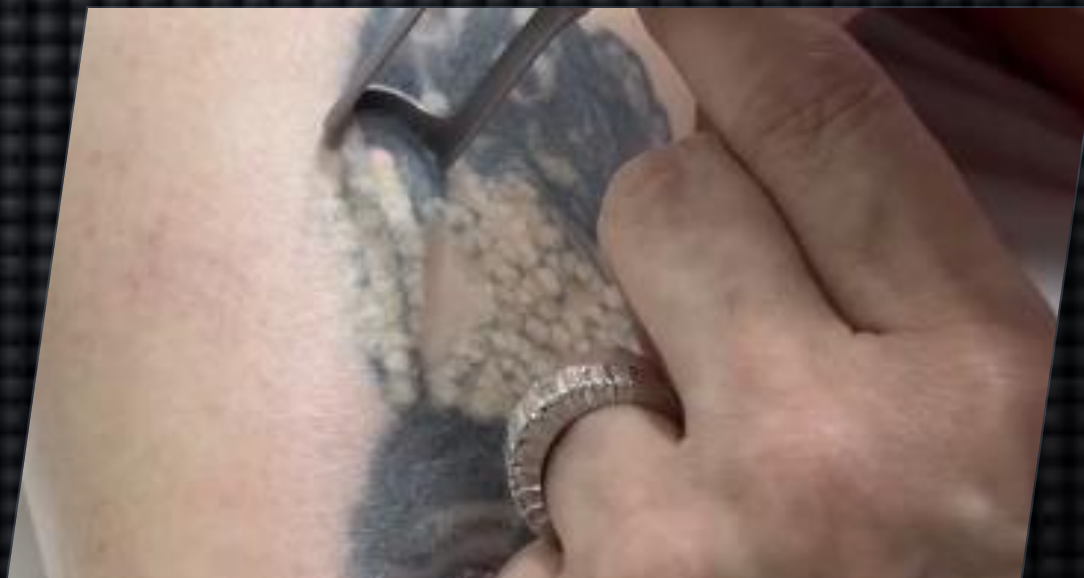


1064nm:
For Dark Color Tattoo And
Pigmentation Removal.

532nm:
For Colorful Tattoo And
Pigmentation Removal




755nm probe

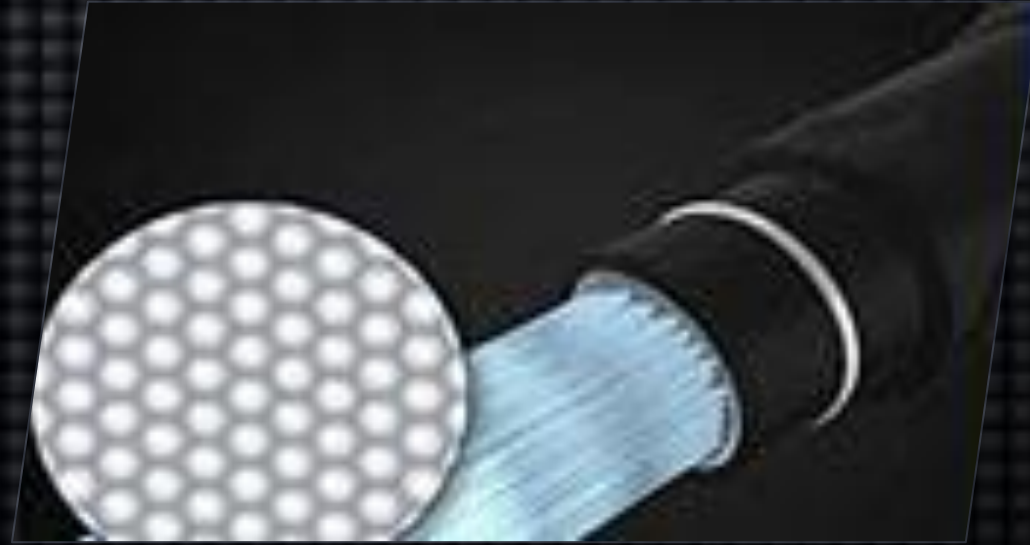

755nm:
Honeycomb Shape,Break The
Melanin,More Dispersed!A Revolution
In Skin Revitalization!




Optional probe

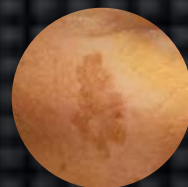

Square Picosecond
Lattice Probe
Suitable for deep
pigment removal, darker
spots, chloasma



Round Honeycomb Probe
It can be used for skin
rejuvenation, skin whitening, rough pore
improvement, aged horniness removal
and blackheads removal.

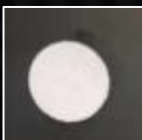


Big Spot Size Probe
Suitable for sweeping, chloasma,
sweeping spots,
sweeping over a large area.

Inside laser



USE CANADA BEAMTECH LASER

- 1) Maintain excellent stability at different energy levels
- 2) Multistage amplification design
- 3) 100w-2000W High power and cooling efficiency

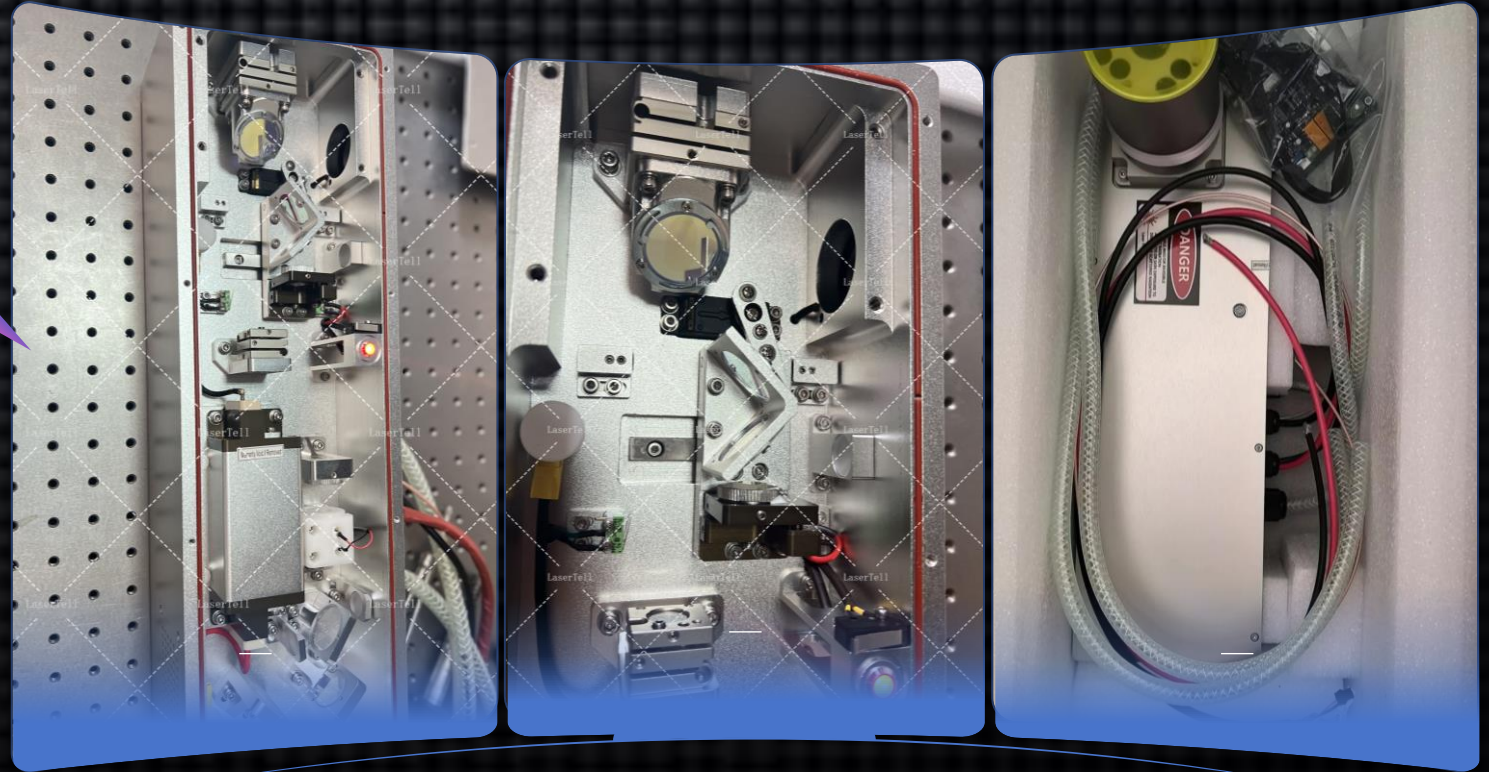




**12.1 inch
color touch screen**

Customized service



Working theory

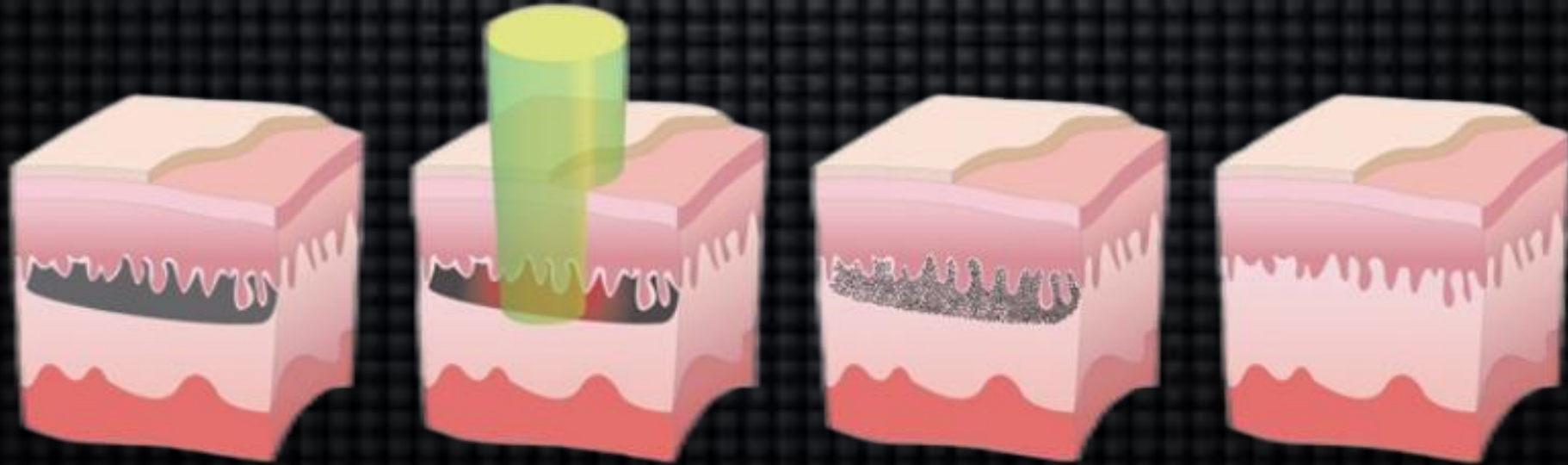


Picosecond laser uses very short pulse output mode, instead of thermal effect, by the principle of light mechanical shock wave, pigment is "shattered" into finely granular via focused energy, are more likely to be absorbed by the body metabolism. Picosecond laser will be to minimize the side effects of thermal effect, can achieve the goal of almost solve all kinds of pigment spots, is better than traditional laser spot whitening effect

Function

Pores




Spot

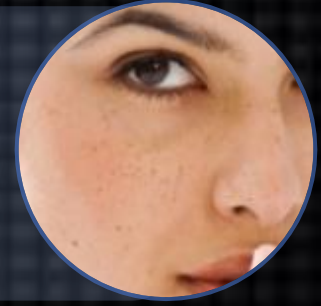

Fine line



Porked scar

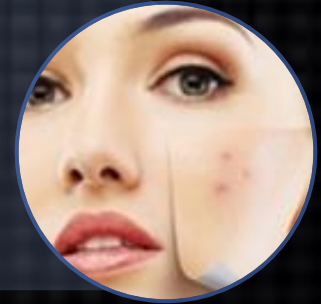


Face dumb



Tattoo



Before&After

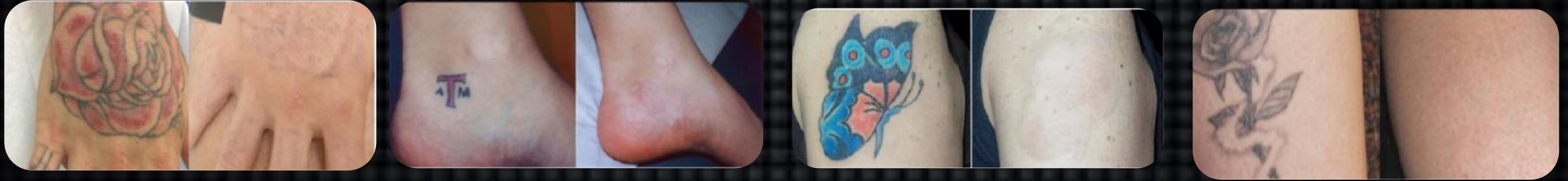

Technical specification

System	Picosecond Q-Switched Nd: YAG Laser
Wave Length	1064nm, 532nm, 755nm
Width of pulse	750ps
Output Energy	1064nm: 100-1200mj; 532nm: 50-600mj
PTP mode	50%energy increasing
Repetition Rate	1-10hz adjustable
Spot Size	φ2-10mm adjustable
Cooling System	Closed-off water circulation +wind cooling
Light Guiding System	Imported 7 jointed nd yag laser articulated arm
Electrica	100-240V, 50/60Hz
Product name	HydroCool_Pico

Real machine display



Ready to ship





Well-recognized by Customers around 80+ Countries abroad with LaserTell Devices



Aftersales service

2 years warranty available (2 +1 charge) ex factory under normal and proper operation; free labor and material but without consumables parts.

Why 2+1 years warranty? (7 Simple Reasons)

1. Normally almost all the medical and aesthetic lasers with only 5 years lifespan, while the first 3 years' guarantee and full maintenance are right the key points.

LaserTell offers 2 years standard warranty and 1 year extension warranty. (Free-2, charge-1)

2. Lifelong Service (no matter 3 years or 5 years, or more...)

3. LaserTell, similar to a local service provider, has a full complement of certified technicians and engineers to assist you with all of your repair and maintenance needs, offering you faster service at much lower prices than the manufacturers

4. 24 hours + 7 days after-sale service and technical supports

5. Service Manual and Safety Instructions available

6. There is a comprehensive training last for 5-7 days in our factory and clinics, including how to operate maintain and adjust the treatment parameters for different treatments, etc