


HYDROCOOL PICOMAX

Product introduction





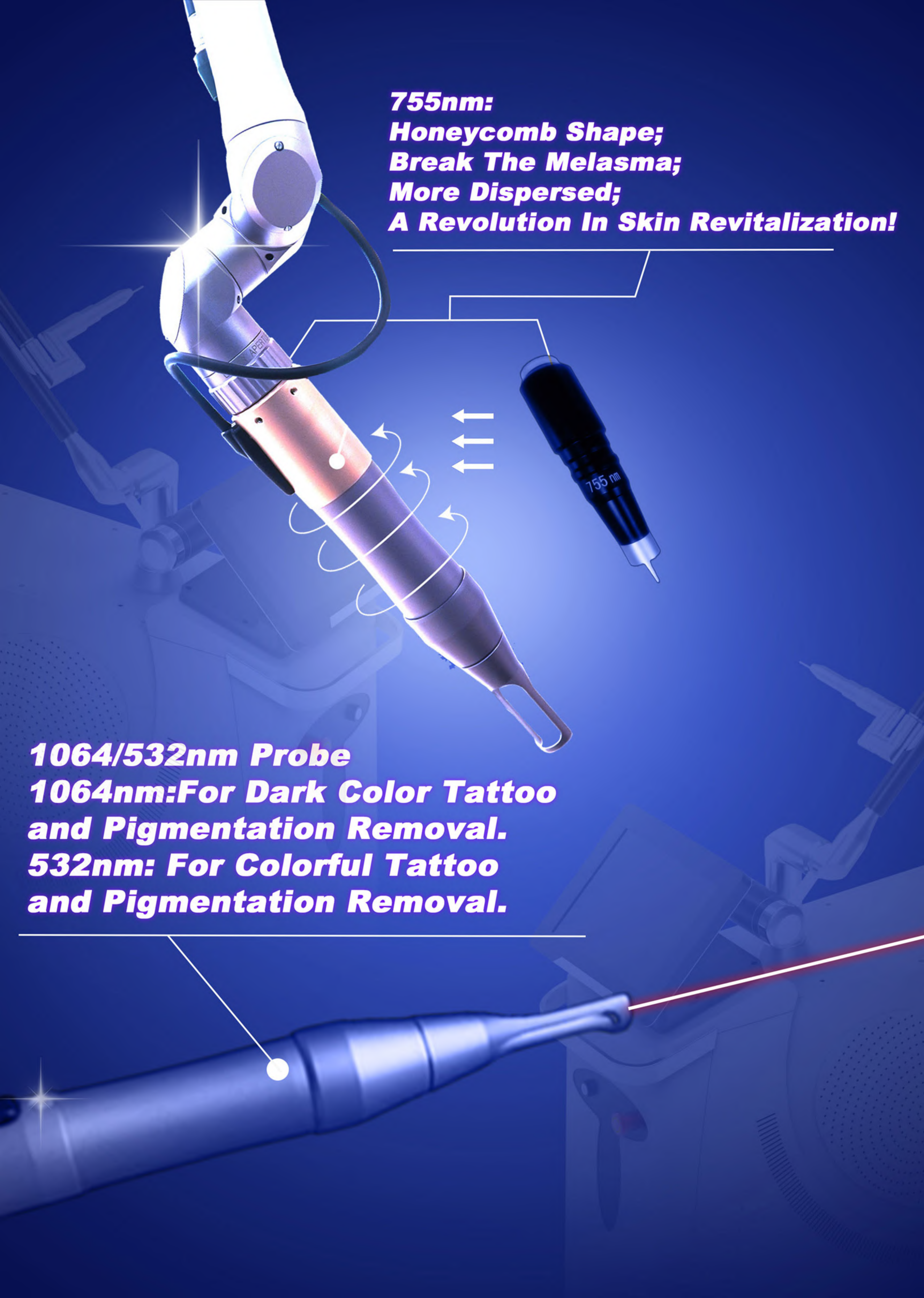


**755nm:
Honeycomb Shape;
Break The Melasma;
More Dispersed;
A Revolution In Skin Revitalization!**

**1064/532nm Probe
1064nm:For Dark Color Tattoo
and Pigmentation Removal.
532nm: For Colorful Tattoo
and Pigmentation Removal.**


Square Lattice Probe
Suitable for deep pigment removal,
darker spots, chloasma.

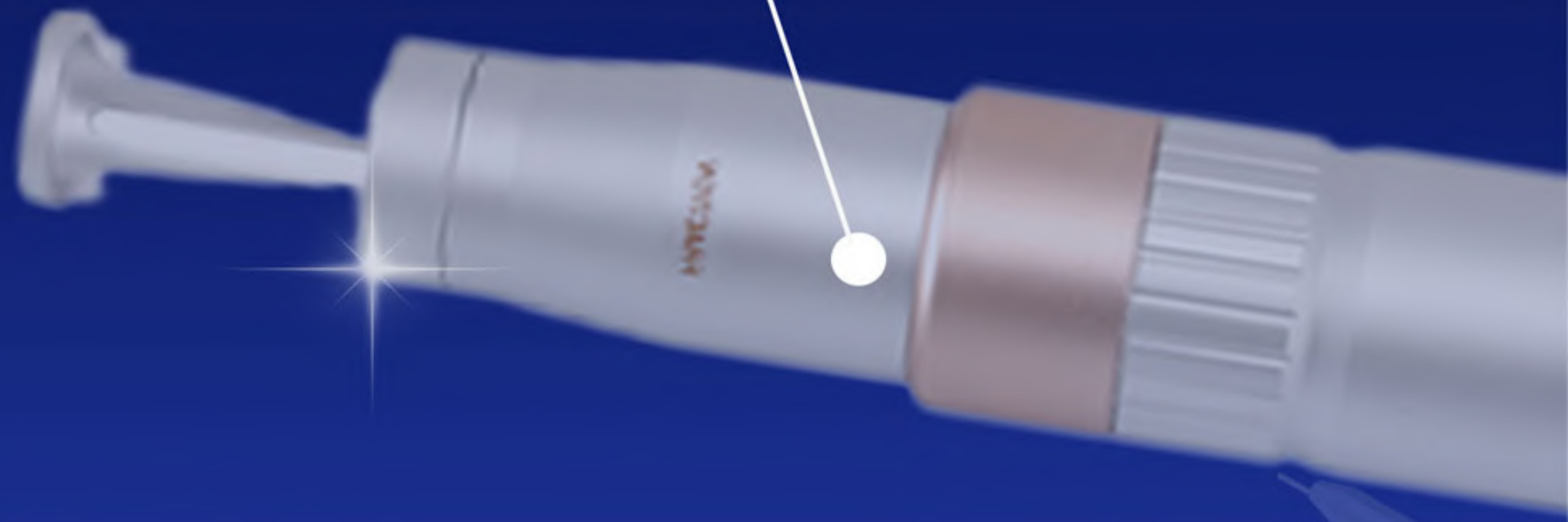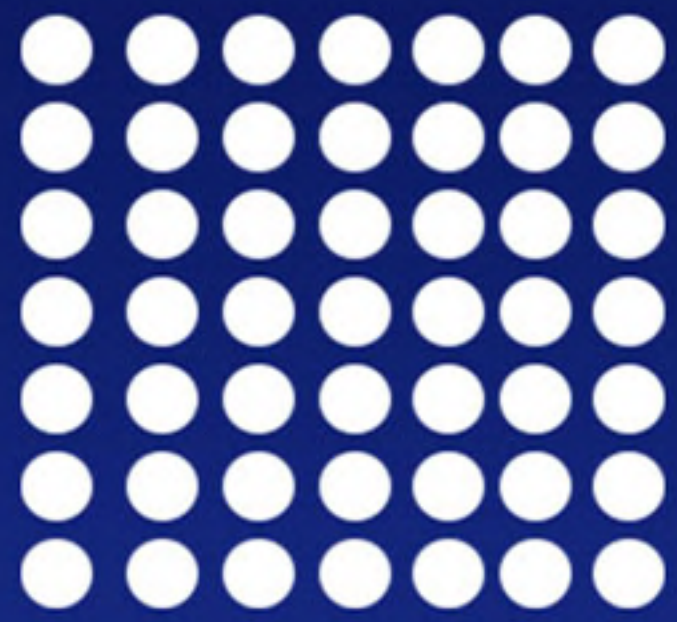

Round Honeycomb Probe
It can be used for skin rejuvenation, skin
whitening, rough pore improvement, aged
horniness removal and blackheads removal.

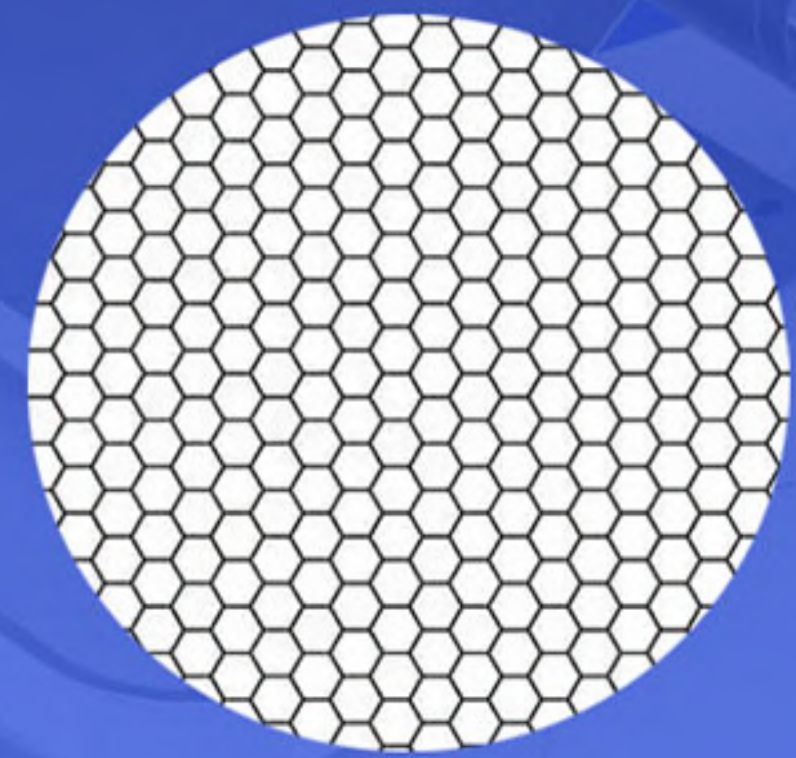

Big Spot Size Probe
Suitable for sweeping chloasma,
sweeping spots, etc.

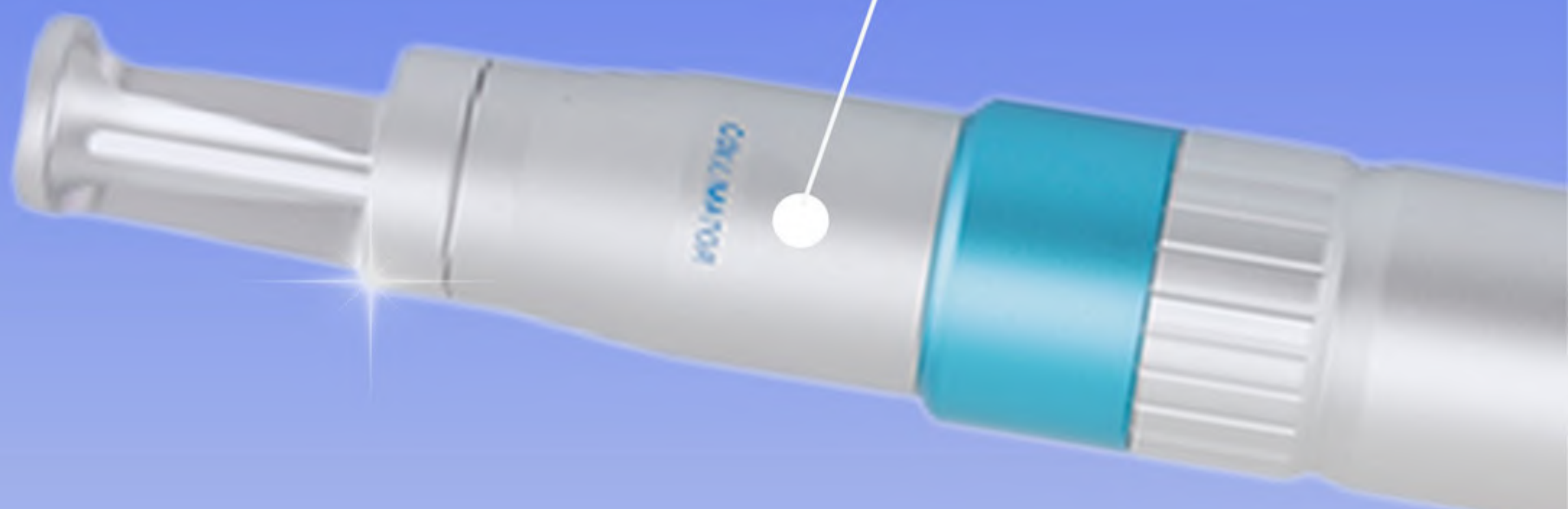

Excellent KOREAN 7

Articular-ARM

to ensure the accuracy of energy output.



Weight balance hammer
keep the arm always be safe

Spot Size: 2-10mm adjustable



Inside placed lens mirror
for better result

FIVE-MODE OPERATION INTERFACE

15.6 inch True-Color Touch Screen

5 modes meet all needs Simple and user-friendly Intelligent system

The latest magical and dynamic interface



INTERFACE INTRODUCTION

532 interface

This function in the process of treatment will emit ultra-long light energy, can quickly lock the skin of the pigment cells, cells will absorb energy after self-exlosion and decomposition, thus improving the skin color, can effectively remove red, green pigment

755 interface

This function can use honeycomb therapy head, can destroy abnormal pigment cells precisely, so that it is gradually cleared by the body, thus achieving the effect of whitening

1064 interface

This function can promote the skin's metabolism, accelerate the old waste of cutin shedding, but also can promote the skin of melanin discharge from the body, help skin more smooth and delicate

Long Pulse mode interface

This function can use black-faced dolls, white porcelain dolls, with skin spots removal, pores shrink, skin texture improvement, skin color uniformity,

PTP mode interface

This function acte in the subcutaneous dermis layer, can promote the reorganization of collagen and elas ticfiders, achieve compact skin, reduce pores, improve the face of small wrinkles and other effects



APPLICATIONS



Pores



Fine lines




Freckles



Spots



Acne scars

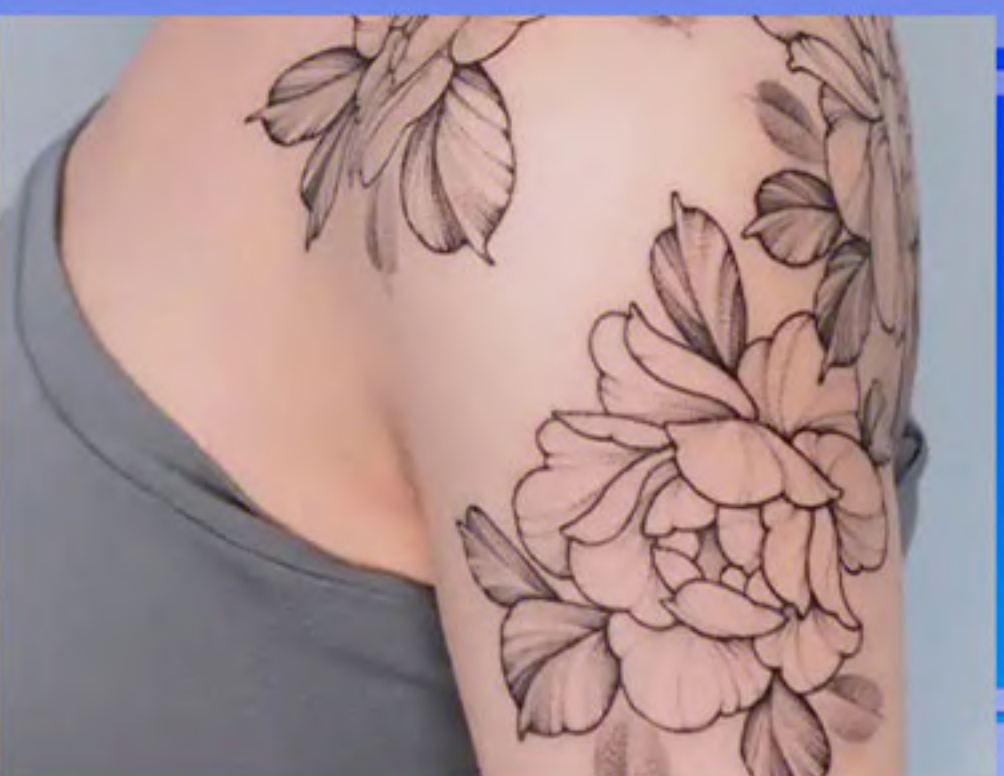


Tattoos

WORKING THEORY

Picosecond Laser uses very short pulse output mode, instead of thermal effect, by the principle of light mechanical shock wave, pigment is "shattered" into finely granular via focused energy, are more likely to be absorbed by the body metabolism.

Picosecond laser will be to minimize the side effects of thermal effect, can achieve the goal of almost solve all kinds of pigment spots.

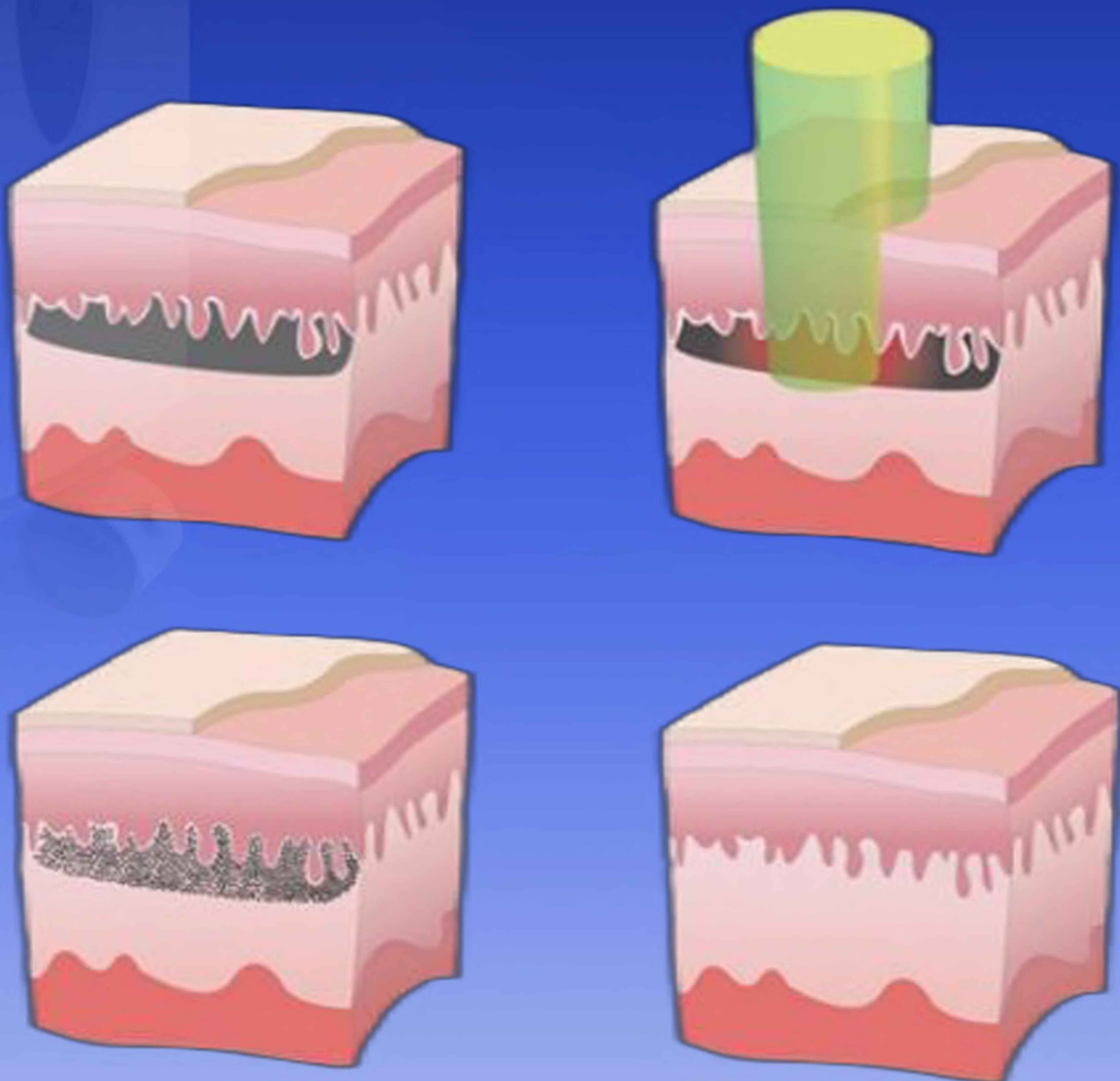

TECHNICAL SPECIFICATIONS

System	Picosecond Laser
Wavelength	1064nm, 532nm,755nm
Pulse Duration	450ps
Output Energy	1064nm:100-1200mJ
	532nm:100-600mJ
PTP & Genesis Mode	100-2400mJ
Repetition Rate	1-10Hz
Spot Size	2-10mm adjustable
Cooling System	Closed-off water circulation + air cooling
Light Guiding System	Korea-imported 7 jointed articulated laser arm
Electrical Requirement	100-240V,50/60Hz
Display	15.6” 4K Touch Screen



PROFESSIONAL CERTIFICATES

AFTER-SALE SERVICE

2 years warranty available (2 +1 charge) ex factory under normal and proper operation; free labor and material but without consumables parts.

Why 2+1 years warranty? (7 Simple Reasons)

1. Normally almost all the medical and aesthetic lasers with only 5 years lifespan, while the first 3 years' guarantee and full maintenance are right the key points.

LaserTell offers 2 years standard warranty and 1 year extension warranty. (Free-2, charge-1)

2 . Lifelong Service (no matter 3 years or 5 years, or more ...)

3.LaserTell, similar to a local service provider, has a full complement of certified technicians and engineers to assist you with all of your repair and maintenance needs, offering you faster service at much lower prices than the manufacturers

4.24 hours + 7 days after-sale service and technical supports

5.Service Manual and Safety Instructions available

6 . There is a comprehensive training last for 5 - 7 days in our factory and clinics, including how to operate maintain and adjust the treatment parameters for different treatments, etc