


iMED_Android

755+808+1064nm Trio

Diode Laser

Android

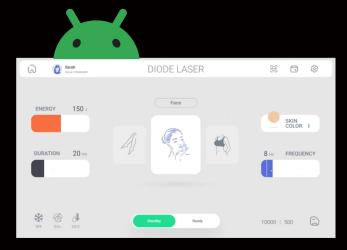
















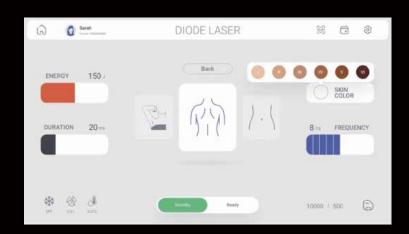






Android system

Diode laser 13.3-inch 4K Android, wifi, Bluetooth, 3D animation smart ipad screen, you can even Google, watch TV, listen to music, and take screenshots.

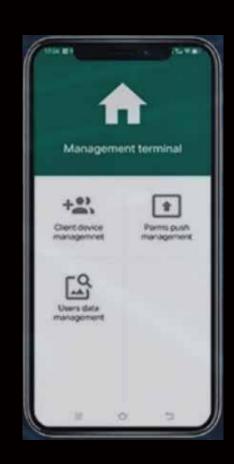

▲ Remote control system▲ Clinical patient database▲ Leasing business system





air instantly once opened:turn off once closed







Remote locking of the machine

Automatic locking at the end of lease time



Patient database Viwe backgeound date remotely







CUSTOMIZATION IN 16 LANGUAGES

Built-in multiple languages to make the operation smoother more languages are being added...















The most Advanced System

Epoch-making tech innovation













Corectell & Sharing better Medicare

10+ Advantage Newly upgraded intelligent system






Start Rental Business

Accrding to the lease agreement set the corresponding lease days and lighting times



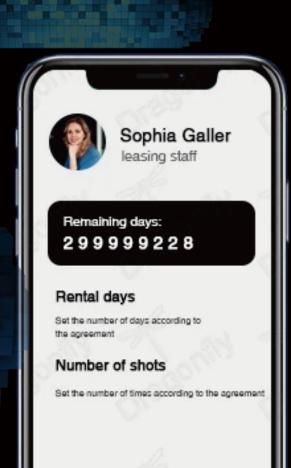




WORKING PRINCIPLE >

DL-808D diode laser hair removal comes with 808nm wave laserstack or 755+808+1064 tripe wave laser stack with semicon.ductors technology





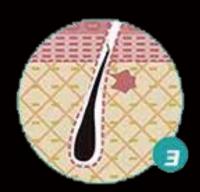



Applicator will emit a laser beam


The laser beam is absorbedby hair follicle melanin.

After absorption, the laser energy will be converted into heat energy. Hair follicles and surrounding stem cells lose their activity and all off So reach hair removale ffect.



What is the Trio wavelength diode laser?

Trio wavelength diode laser proprietary hair removal innovation is clinically proven method of laserhair removal that is virtually painless, and it is the safe laser hair removal methodfor all skin types, including dark and tanned skin. Its unique gradual heating andrapid/fast pulse delivery method effectively damages the hair follicle and hindersre-growth, while preventing injury to the surrounding tissue.

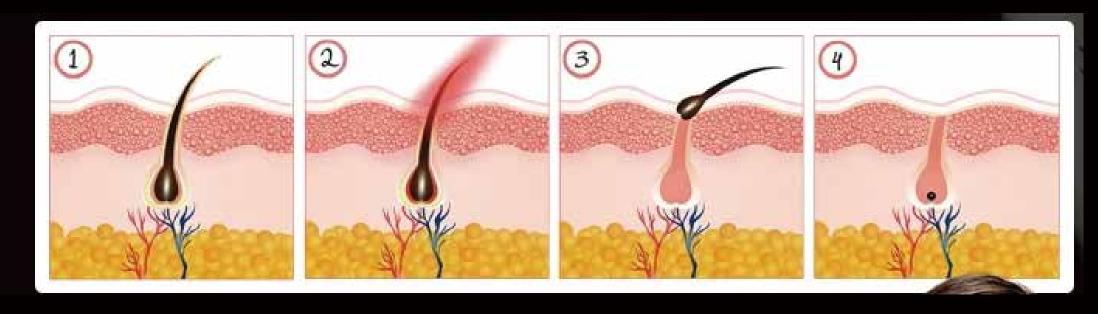


Theory



The principle of diode laser hair removal treatment is based on selective photothermal decomposition. The melanosomes in the hair follicle can selectively absorb the energy of the laser. Thelaser energy emitted by the machine is easily absorbed by the colored hair follicle without damaging the epidermal tissue. The energy emitted by the laser is absorbed by the hair and the pigment in

It is competed by the state thus increasing the temperature of the hair follicle. When the temperature rises to a certain level, the hair follicle will be irreversibly damaged, and the hair will lose its original environment and be completely removed.

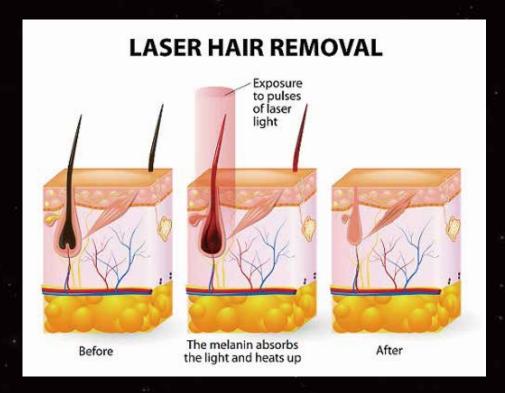





Coherent laser bar Import from USA 50 million shots!

More laser shots-50 millions Real 3 wavelength 755 +808+1064 High power each bar 100w Stable quality longer lifetime












Handle

Multifunctional Detachable Design

Before & After





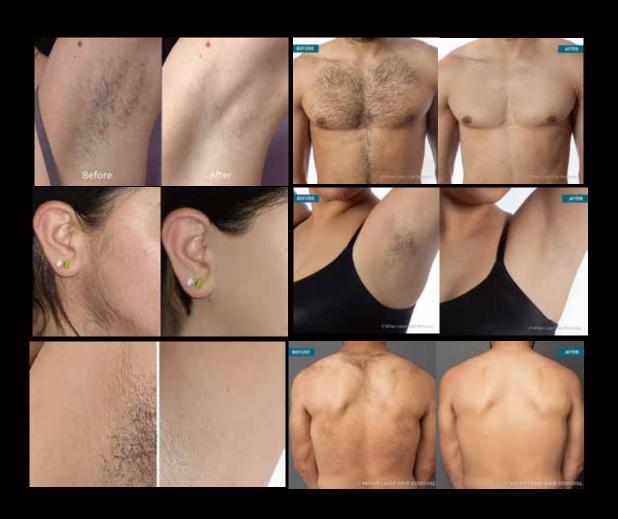





Model Name

Machine Power

Laser Stack Power

Diode EnergyDensity

Spot Size

Frequency

Display

Net Weight

Dimensions

Stand-by Working

iMED Andriod

3000W

2000W

1-150J/cm²

1-10Hz

Electrical Requirements 110/230VAC,15/20A max,

13.3" HD Android screen

28KG

30*50*40cm

Continuously 24 hours



Shoring better Meditore

Views of TruMEDiO7















Shadon better Medicare

Our Certificates























Well-recognized by Customers around 80+ Countries abroad with LaserTell Devices





After-sale Service

2 years warranty available (2 + 1 charge) ex factory under normal and proper operation; free labor and material but without consumables parts.

Why 2+1 years warranty? (7 Simple Reasons)

- 1. Normally almost all the medical and aesthetic lasers with only 5 years lifespan, while the first 3 years' guarantee and full maintenance are right the key points.
- LaserTell offers 2 years standard warranty and 1 year extension warranty. (Free-2, charge-1)
- 2. Lifelong Service (no matter 3 years or 5 years, or more...)
- 3. LaserTell, similar to a local service provider, has a full complement of certified technicians and engineers to assist you with all of your repair and maintenance needs, offering you faster service at much lower prices than the manufacturers
- 4. 24 hours + 7 days after-sale service and technical supports
- 5. Service Manual and Safety Instructions available
- 6. There is a comprehensive training last for 5-7 days in our factory and clinics, including how to operate maintain and adjust the treatment parameters for different treatments, etc